

WANJIE PULSE LAVAGE SYSTEM

A fully disposable instrument,designed for use in orthopaedic surgical procedures,trauma and general wound care.

The WANJIE pluse is a versatile instument which combines convenience and power for removal of necrotic tissue,bacteria and foreing material whilst providing controlled pressure for efficient cleansing and in a double sterile pack and includes a lightweight hand piece,two tips,one of which is a dual cone splash shield tip,designed for optimal fluid containment,therefore helping to control the dangers associated with transient spray.



Key Features Lightweight Design	Sleek,ergonomic design comfortable to handle and use.
Smooth Action Toggle Trigger	Fingertip control for either high or low speed fluid delivery.
Secure Tip Locking Mechanism	Ensures secure attachment of the tips.
Textured Barrel and Handle	Provides a secure and comfortable grip, ensuing greater control operative procedures.
Independent Fluid and Suction Tubes	Provides concentric fluid delivery and suction.



BEIJING WANJIE MEDICAL DEVICE CO.,LTD.


Ad: Building 8-2,Zhongguancun,M-park, Yongwang West rd,
Daxing District,Beijing, China
Tel: 010-56358538, +86 13801320400
FAX: 010-88455330
E-mail: zhangmei@wanjiety.com.cn
Web: www.wanjiety.com.cn



WANJIE Pulse Lavage System



BEIJING WANJIE
MEDICAL DEVICE CO.,LTD.

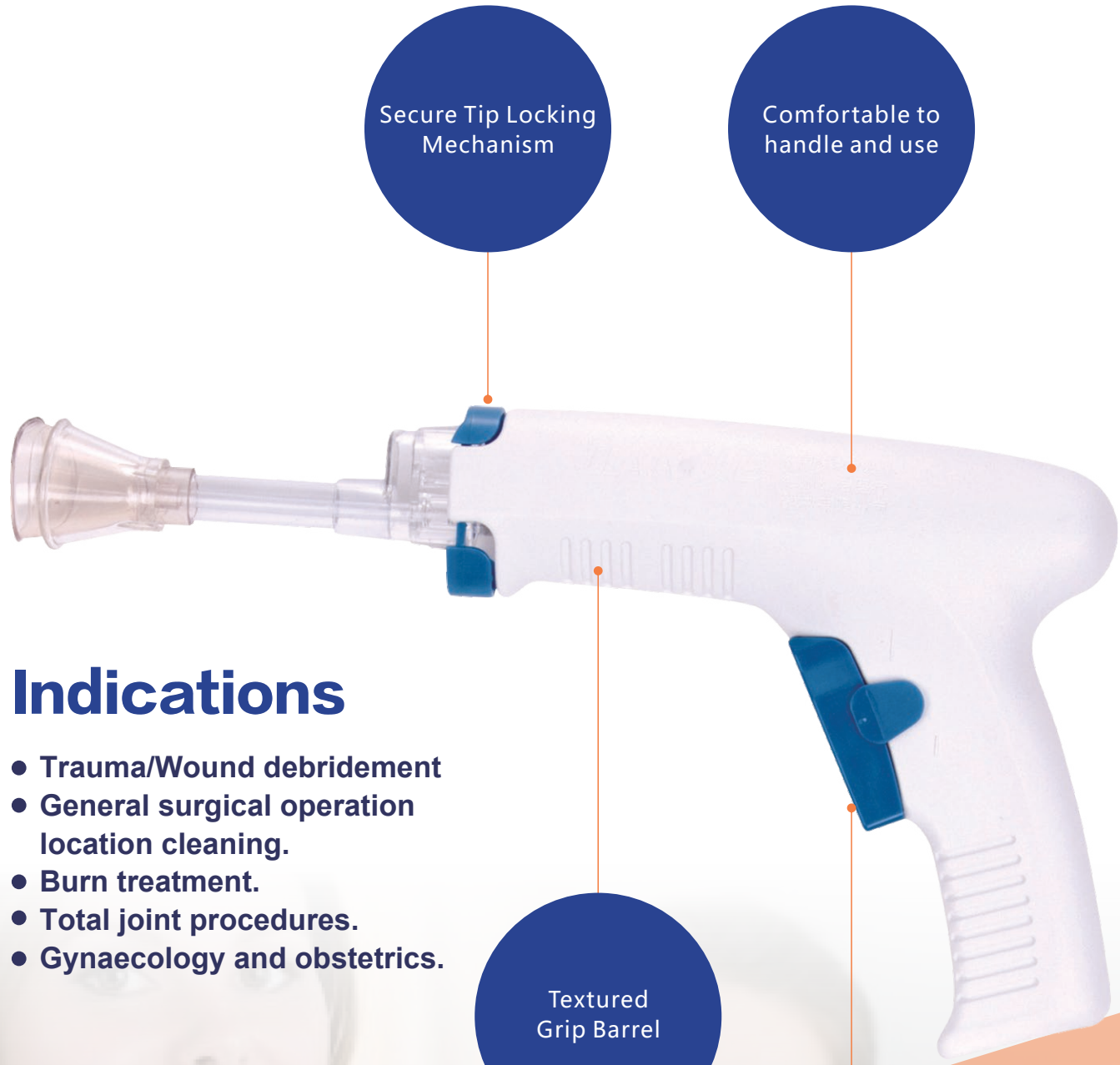

Indications

- Trauma/Wound debridement
- General surgical operation location cleaning.
- Burn treatment.
- Total joint procedures.
- Gynaecology and obstetrics.

Toggle Trigger for high or low speed fluid delivery

WANJIE Pulse Lavage System


Pulse Lavage Drainage Device AC

The Pulse Lavage Drainage Device AC is a disposable, environmentally friendly system that includes a handle gun, a high capacity fan spray debridement tip, a high capacity narrow intra modularly tip, and a reusable AC power pack that eliminates the need for batteries. Our Pulse Lavage Drainage Device AC saves your facility money by eliminating the expense of battery waste all while providing a better debridement and wound care experience. It provides versatility and reliability for effective, thorough and environmentally friendly debridement.



What less could mean

A hospital that utilizes 1,000 Pulse Lavage Drainage Device units each year will pay for the disposal of approximately 600 pounds of battery waste. Over a five year period, this is the equivalent to throwing away a full size car weighing 3,000 pounds (1.5 tons).

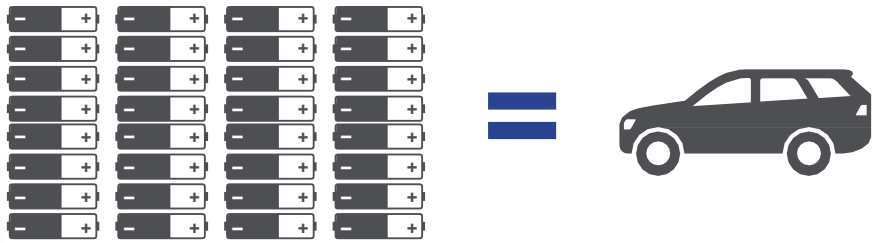

Pulse Lavage Drainage Device

The pulse lavage drainage device is used for wound debridement. It combines convenience and power for a smooth and efficient flow of pulsatile irrigation and debridement. The major clinical department is orthopedic department. It can be used for hip and knee replacement, trauma debridement etc. The Pulse Lavage Drainage Device operates on eight AA batteries.

Tip Guide

Name	Tip	Tip Description	Spray Pattern	Feature
Narrow Intra Modularly Tip		High Capacity Single Stream Tip w/ Splash Shield, 11.9 cm length		"Pinpoint" Spray Single FluidStream
Fan Spray Debridement Tip		High Capacity Fan Spray Tip w/ Splash Shield, 11.9 cm length		Wide Spray Pattern
Splash Shield		Dual-cone		Designed to minimize splash back and help control the dangers associated with transient spray
Battery Pack		Type BF applied part, contains eight AA batteries		Easy to carry and without electric shock
Reusable AC Adapter		Comply with EU standards, voltage: AC 220V		Energy saving, environmental protection, low consumption and reusable

Ordering Information

Pulse Lavage Drainage Device AC	WJPL0006001	Handle Gun
WJ-06 AC	WJPL0000001	High Capacity Fan Spray Debridement Tip
	WJPL0000002	High Capacity Narrow Intra Modularly Tip
	WJPL0006002	Reusable AC Adapter
Pulse Lavage Drainage Device	WJPL0006101	Handle Gun
WJ-06	WJPL0000001	High Capacity Fan Spray Debridement Tip
	WJPL0000002	High Capacity Narrow Intra Modularly Tip
	WJPL0006101	Battery Pack
Pulse Lavage Drainage Device	WJPL0017001	Handle Gun
WJ-17	WJPL0000001	High Capacity Fan Spray Debridement Tip
	WJPL0017001	Battery Pack
Pulse Lavage Drainage Device	WJPL0018001	Handle Gun
WJ-18	WJPL0000001	High Capacity Fan Spray Debridement Tip
	WJPL0000003	Splash Shield
	WJPL0018001	Battery Pack

Contact your Beijing WANJIE representative or visit us at www.wanjiety.com.cn

This documentation is intended exclusively for physicians and is not intended for laypersons. Information on the products and procedures contained in this document is of a general nature and does not represent and does not constitute medical advice or recommendations. Because this information does not purport to constitute any diagnostic or therapeutic statement with regard to any individual medical case, each patient must be examined and advised individually, and this document does not replace the need for such examination and/or advice in whole or in part. Please refer to the package inserts for important product information, including, but not limited to, indications, contraindications, warning, precautions, and adverse effects.